

## Occlutech<sup>®</sup> Delivery Set III (ODS III)

Intended for the introduction of diagnostic  
and therapeutic devices into the vasculature

Designed for Optimal Delivery

# Occlutech<sup>®</sup> Delivery Set III

ODS III


The Occlutech Delivery Set III is designed for optimal performance while providing exceptional pushability, trackability and kink resistance.

The Occlutech Device consists of a braided sheath with a three-way stopcock side port, hemostatic valve, dilator and loader. The distal tip of the sheath is radiopaque and contains two side flush holes. The dilator can be locked onto the hub of the sheath.


**A** Dilator

**B** Sheath

**C** Loader

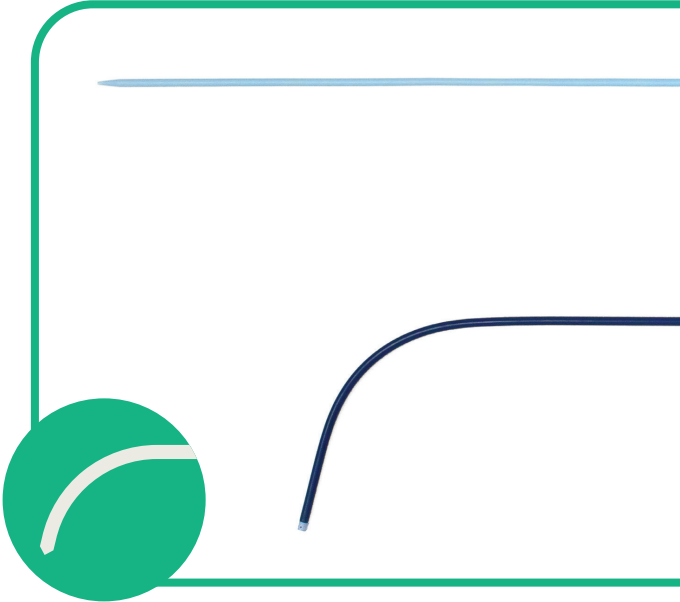



45°

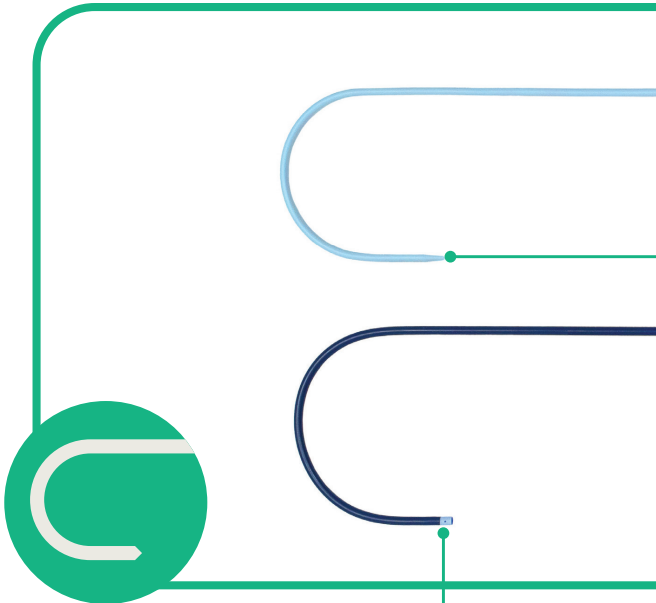

180°



Distal marker band and side holes designed for accurate positioning and flushing

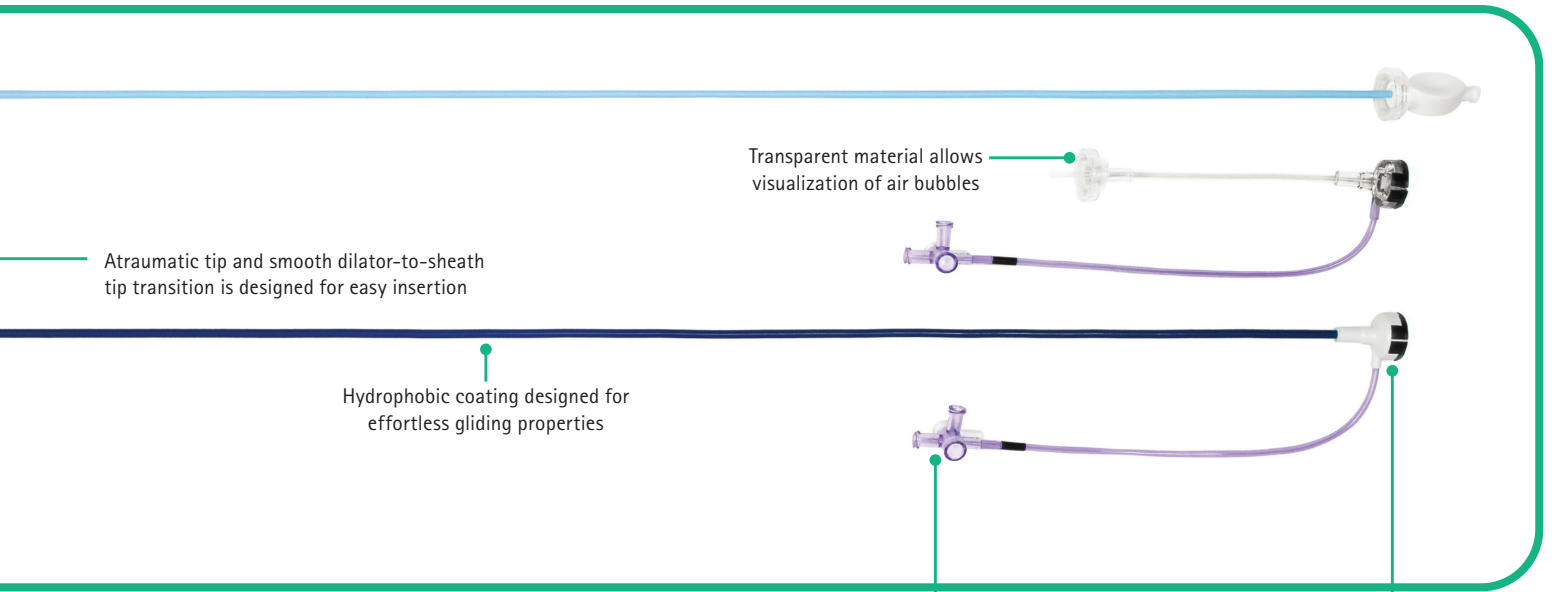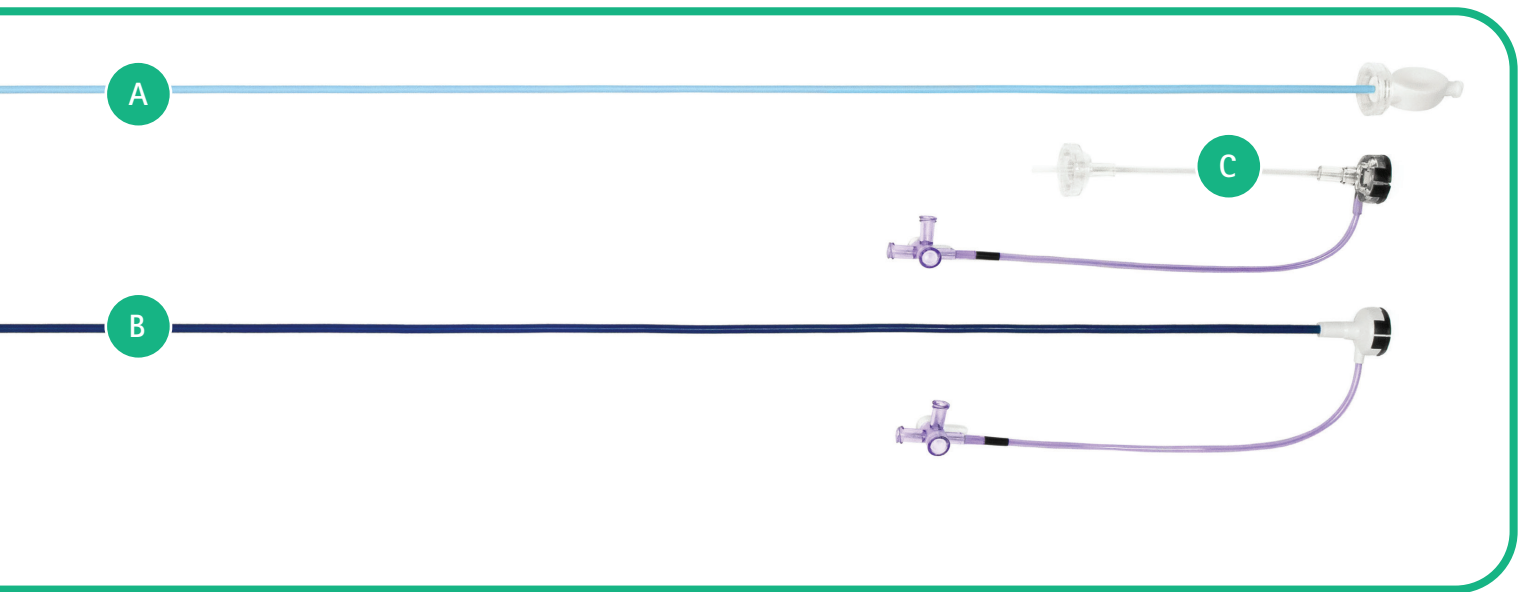






















Three-way stopcock



Hub with hemostatic valve designed to minimize blood loss

# Occlutech® Delivery Set III

## Ordering Information

| Sheath Type   | Reference Number   | Size (Fr) | Usable Length Sheath (mm) | Usable Length Loader (mm) | Distal Angle (°) | French Size Color Code  |
|---|--|-----------|---------------------------|---------------------------|------------------|---|
| 45°<br>Regular<br>Length<br><br> | 98DS006*   | 6         | 790                       | 85                        | 45               |    |
|   | 98DS007  | 7         | 790                       | 85                        | 45               |    |
|   | 98DS008*   | 8         | 790                       | 125                       | 45               |    |
|   | 98DS009  | 9         | 790                       | 125                       | 45               |    |
|   | 98DS010  | 10        | 790                       | 125                       | 45               |    |
|   | 98DS011  | 11        | 790                       | 165                       | 45               |    |
|   | 98DS012  | 12        | 790                       | 165                       | 45               |    |
|   | 98DS014*   | 14        | 790                       | 165                       | 45               |    |
| 45°<br>Long<br>Length<br><br>  | 98DL006*   | 6         | 1100                      | 85                        | 45               |    |
|   | 98DL007*   | 7         | 1100                      | 85                        | 45               |  |
|   | 98DL008*   | 8         | 1100                      | 125                       | 45               |  |
|   | 98DL009*   | 9         | 1100                      | 125                       | 45               |  |
|   | 98DL010*   | 10        | 1100                      | 125                       | 45               |  |
|   | 180°<br>Regular<br>Length<br><br> | 98US006*  | 6                         | 790                       | 85               | 180   |
| 98US007*  |  | 7         | 790                       | 85                        | 180              |  |
| 98US008*  |  | 8         | 790                       | 125                       | 180              |  |
| 98US009*  |  | 9         | 790                       | 125                       | 180              |  |
| 98US010*  |  | 10        | 790                       | 125                       | 180              |  |
| 98US011*  |  | 11        | 790                       | 165                       | 180              |  |

\* Made to order

The Occlutech Delivery Set III is intended for the introduction of diagnostic and therapeutic devices into the human vasculature, including but not limited to intracardiac, renal or other placements.

Refer to the Occlutech Delivery Set III Instructions for Use for complete prescribing indications, contraindications, warnings and precautions. Occlutech is a registered trademark of Occlutech Holding AG.

For more information, please contact your Medical Device Specialist or call 877-836-2228.

Rx only.

Distributed by:

B. Braun Interventional Systems Inc. | a member of the B. Braun Group of Companies

Bethlehem, PA 18018 | USA | 877-836-2228 | www.bisusa.com

Technical Support: 800-443-8362



Not made with natural rubber latex



Not made with DEHP



Not made with PVC